



AccessRx Exchange

Helping to Close the Loop on Medication Administration.

The AccessRx Series of medication delivery carts provides solutions to more securely and efficiently transport and administer medications at the patients' bedside.



Maximum Number of Bins	12	24	12
Configurable Drawer Options	\checkmark	\checkmark	\checkmark
Exchangable Medication Cassettes		\checkmark	
Individual Locking Patient Drawers	\checkmark		\checkmark
Software links to the HIS	\checkmark		
Controlled Substance Drawer Option	\checkmark	\checkmark	\checkmark
Patient Bin Assignment	Automatic	Manual	Manual
Antimicrobial Infused Cart Material		*does not include cassette	

KEY BENEFITS



MEDICATION DELIVERY WORKSTATIONS

Closing the loop of medication delivery between pharmacy and the patient bedside.

AccessRx MD








Medication Security & Visibility at the bedside

- Individually locking, patient-specific drawers can be accessed by scanning the patient ID barcode
- Real-time bedside medication administration, barcode patient identification and direct eMAR API access streamlines BCMA process, reduces redundancy and increases accuracy
- Patient drawers automatically assigned using the ADT (Admission/ Discharge/Transfer) interface, eliminating time consuming and error-prone manual drawer labeling
- Utilizes the hospital's Active Directory for consistent user authentication
- Configurable drawer display fields allow customization of the user interface

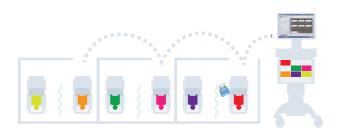

Streamline Your Workflow

- Secures and tracks patient specific medications for up to 12 patients in individually-locking drawers, reducing trips and time spent at the ADC
- ADT and patient orders transmitted via wireless network to the cart to facilitate faster medication delivery
- Reduce trips to and time spent at the automated dispensing cabinet



More Control & Flexibility for Pharmacists across a variety of pharmacy distribution models

- Browser-based AccessManager[™] console provides access to reports while allowing for Active Directory-compatible user management, fleet management and system configuration
- Configurable, time-based access controls reduce the risk of medication administration at an incorrect time and aid in pharmacy reconciliation for unused medications
- Optional medications orders interface Rx MD can act as a patient medication order filler and verify that an order has been administered.
- Secure Pharmacy email alerts for transferred & discharged patients



Scalable & Centralized Software Architecture and Integration

- Multi-Org scaleable architecture ability to scale from two bed clinic to multi facility health system
- Ability to set up health system access rights, such as by department
- Enterprise reporting available
- Enterprise multi-facility deployment
- HL7 compliant, ADT interface (patient demographics), and RDE interface (Pharmacy approved orders)
- Configurable interface engine easily communicates with the HIS
- Project management, installation and training by TouchPoint Medical's professional services team

Configure your system to meet your medication and supply requirements.

		102mm (4in) Drawer Specifications	152mm (6in) Drawer Specifications	Controlled Substance Full Width Drawer Specifications
	Dimensions Width Height Depth Volume	91.4mm (3.6in) 238.8mm (9.4in) 55.9mm (2.2in) 1213cm³ (74in³)	142.2mm (5.6in) 238.8mm (9.4in) 55.9mm (2.2in) 1885cm ³ (115in ³)	271.8mm (10.7in) 228.6mm (9in) 50.8mm (2in) 3195cm³ (195in³)

MEDICATION DELIVERY WORKSTATIONS

Driving efficiency & security with your centralized medication distribution.

AccessRx Exchange











AccessRx Exchange[™] Workstation on Wheels

- Restock easily wherever your need. Can easily be filled in the pharmacy and exchanged or directly at the cart.
- Key lock or eLock security options
- Easily retrofit your existing AccessPoint[™] workstation on-site with the Exchange Cassette system
- May be added to powered or nonpowered carts
- Works with both manual and electric lift carts


Exchangable Cassettes

- Restock easily and directly on the mobile workstation or wall mount
- Available in two, three or four tier sizes & multiple bin sizes and supply drawer options
- Load rated up to 2.2kg (5.5lbs) per tier
- Compact footprint fits 600mm x400mm cart or (ISO 3394 modular systems) cabinet

Transfer Carts

- All transfer carts include locking tambour door
- All transfer carts include extended handle
- Lock options include key lock or e-lock
- Single width, front loading option
- Double width, side loading option



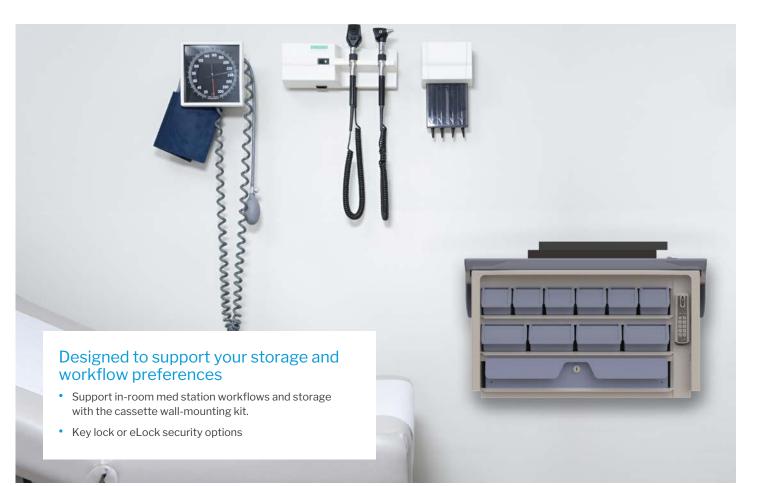





Easily adapt your existing AccessPoint[™] fleet with the AccessRx Exchange[™] system

- Retrofit your existing workstations on-site with minimal downtime
- May be added to powered or nonpowered carts
- Works with both manual and electric lift carts



Flexible Locking Options

Push Button eLock Features:

- 20 codes per lock: 19 User codes and one Supervisor code; code length is 4-8 digits
- 7500 cycles (extended battery life)
- ADA compliant
- Low battery indicator
- Battery compartment is easily accessible from outside of the cassette

Key Lock Features:

- Removable lock core
- Unique lock/key series for cassette carrier, cassettes, and controlled substance drawer for enhanced security
- Side shuttered lock with symmetrical keying

MEDICATION DELIVERY WORKSTATIONS

Store and secure your medications to the bedside.

AccessRx Secure






Secure your medications throughout your workflow

- Individual locking drawers
- Easy-to-use software
- Storage configurations to secure a variety of medications and controlled substances
- Smoked polycarbonate 4" & 6" drawers
- Aluminum locking storage and controlled substance drawers
- Configurable drawer display fields allow customization of the user interface
- Add AccessRx MD[™] software options to existing AccessRx Secure [™]carts for a scalable, custom solution

Configure your system to meet your medication and supply requirements.

With a wide offering of sizes and locking options, you can have peace of mind knowing your medications and supplies are securely stored to your requirements.

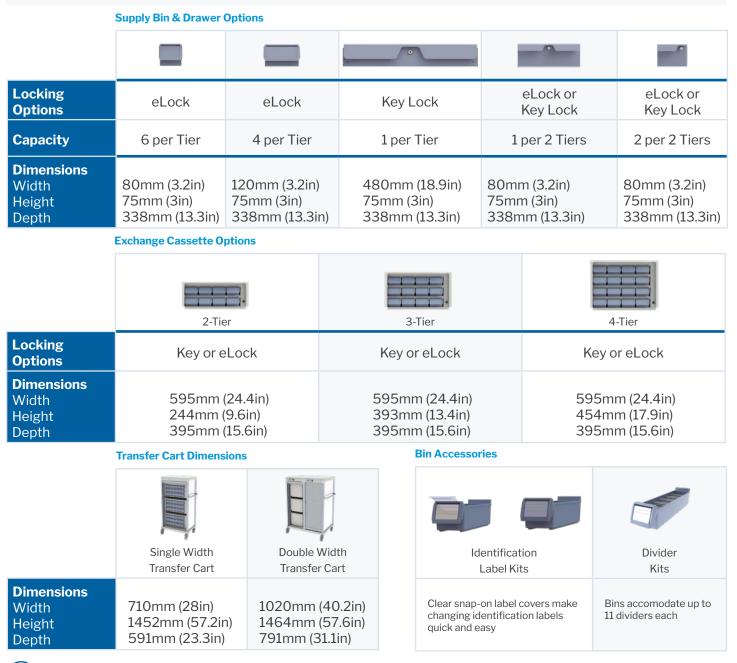

Contact your TouchPoint Medical representative to learn more and discuss additional accessories specific to your point-of-care application. Visit **touchpointmed.com/contact** to get in touch with your local representative.

20200121-BR-MedicationDeliveryCarts-EN

